



Naturopath Hydrogen Generator 150ml/min Portable Inhalation Machine

Our Product Introduction

for more products please visit us on weizi-medical.com


Basic Information

- Place of Origin: Guangdong, China
- Brand Name: Weizi
- Minimum Order Quantity: 1
- Price: \$185-205




Product Specification

- Production Rate: 2000pcs/month
- Weight: 1.65kg
- Dimensions: 154*154*241mm
- Hydrogen Output: 150±10 ML/min
- Oxygen Output: 75±10 ML/min
- Flow Rate: 100, 150, 225 ML/min
- Timer Setting: 1,2,3 Hours
- Working Life: 5,000 Hours
- Power: 65w
- Voltage: AC90-260V, 50/60Hz
- Water Replacement: Every 20 Hours
- Water TDS: <3ppm
- Daily Usage: 1/2 Hour Initially
- Max Usage: 1 Hour Twice Daily
- Technology: SPE And PEM

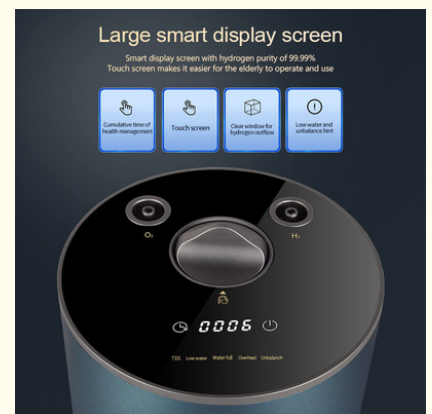

More Images

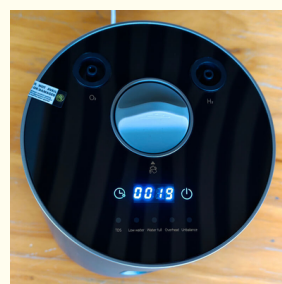

Product Description

150ml Hydrogen Power Inhalation Machine

Portable Molecular Hydrogen Breathing System

Molecular hydrogen is a small, powerful and effective antioxidant which specifically targets harmful free radicals by neutralising them via electrons without any adverse byproducts.

**Portable Oxy Hydrogen Generator
for the loved**

Hydrogen purity up to 99.99%
High-quality hydrogen production

- Antioxidative anti aging
- Enhancing immunity
- Relieving fatigue
- Relaxing mood
- Helping sleep
- Beautifying skin
- Anti inflammatory

Advantages of Our Hydrogen Breathing System

The smallest molecule has its advantages. This advantage allows hydrogen to diffuse into every organ and cell in the body which allows maximum availability for **achieving antioxidant, anti-inflammatory and anti-allergy effects**.

5 detection systems for the loved one



TDS value detection

Ensure that TDS (total dissolved solids in water quality) < 10 ppm to prevent the service life from being reduced due to blockage inside the device



Overflow detection

Keeping the water line normal at all times to prevent overflow



Low water detection

Preventing dry burning and prolonging the service life of the device



Overheat detection

Preventing wear and tear to the device caused by excessive temperature



Unbalance detection

Preventing wear and tear to the device caused by dumping



Recommended Usage Frequency

Hydrogen Inhaler is quiet and straightforward to use. Just add distilled water, plug it in and inhale hydrogen at your convenience.

Use the single inlet cannula connected to the Hydrogen outlet


Start with **30 minutes per day** initially

After a few days, increase to **30 minutes in the morning and 30 minutes in the evening**

After 1-2 weeks, progress to **1 hour in the morning and 1 hour in the evening**

This break between sessions allows time for your body to return to its normal baseline levels.

Hydrogen inhalation+oxygen inhalation
Oxy hydrogen mixture
One device serves several purposes

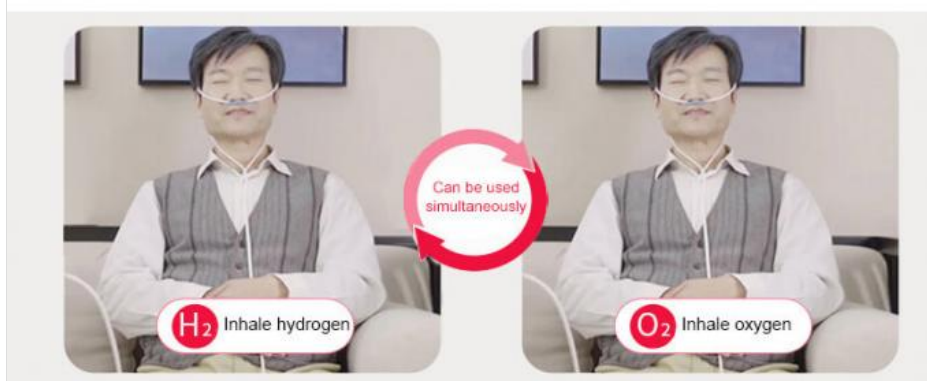

Operation & System Maintenance

Incorporating intelligent sensor systems to deliver our hydrogen breathing system with continual purity analysis, adaptive safety programs and automatic built-in cleaning for everyday ease.

Operating instructions

Operate with ease

The operation is convenient, and users can easily enjoy hydrogen health in a few simple steps



① Connect the power supply



② Remove the water tank cover



③ Add Yibao purified water



④ Remove the sealing plug
Drain for 10-15 seconds



⑤ Connect the suction pipe



⑥ Inhale hydrogen

Technical Specifications



Technology	The latest generation SPE and PEM technology
Flow Rates	100, 150, and 225 ml/min (total flow rate)
Max Hydrogen (H2) Output	150±10 ml/min
Max Oxygen (O2) Output	75±10 ml/min
Timer Settings	1, 2, 3 hours
Working Life	5,000 hours
Maintenance	Replace distilled water every 20 hours of operation
Power	65W
Voltage	AC90-260V, 50/60Hz

Note: Always use high-quality Distilled Water with Total Dissolved Solids (TDS) of less than 3ppm.

Applications

Applications

Safe and comfortable, suitable for all kinds of places



Home



Hotel



Gym



Travelling



Airport



High-speed railway station

About Weizi Medical

Shenzhen Weizi Medical Co., Ltd is a high-tech enterprise, focusing on the R&D, production, and sales of hydrogen-ion health products such as hydrogen-ion inhalation machines, hydrogen-ion water machines, oxygen and hydrogen-ion tanks and etc. Some of our products have got CE, FCC, VOC, UKCA, CNAS and RoHS certificates.

Since 2012, we have collaborated with more than 10 hospitals in China. During the 2019 pandemic, our hydrogen-ion generator machines proved highly effective. Weizi has its own core patented technology and a team of medical professors, giving us strong R&D capability to continuously develop negative hydrogen-ions medical treatment machine devices.

We provide professional ODM/OEM project services and are an anion environment solution provider. With a professional overseas sales team and one of our shareholder doctors long-term working in Germany, our aim is to explore the source of life.



Shenzhen Weizi Medical Co., Ltd.



+8617376733796



lzy@szwzyl.com



weizi-medical.com

2L, 2nd Floor, Building 1, Quanzhi Science and Technology Innovation Park Houtingshe, Shajing Street, Baoan District, Shenzhen, Guangdong Province

